REGLO *Digital*

Schlauch- pumpe

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Schutzgrad
IP 30

2 Kanäle

ISM 831	6 Rollen
ISM 832	8 Rollen
ISM 596	12 Rollen

4 Kanäle

ISM 833	6 Rollen
ISM 834	8 Rollen
ISM 597	12 Rollen

Tubing pump

Microprocessor
controlled

Protection rating
IP 30

2 channels

ISM 831	6 rollers
ISM 832	8 rollers
ISM 596	12 rollers

4 channels

ISM 833	6 rollers
ISM 834	8 rollers
ISM 597	12 rollers

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Classe de protection
IP 30

2 canaux

ISM 831	6 galets
ISM 832	8 galets
ISM 596	12 galets

4 canaux

ISM 833	6 galets
ISM 834	8 galets
ISM 597	12 galets

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Deutsch

New
with menu settings

English

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⚠ Please note

We recommend you to read this operating manual carefully.

When operating a pump, certain hazards cannot be excluded.

ISMATEC SA does not take liability for any damage resulting from the use of an ISMATEC® pump.

ISMATEC SA does not admit responsibility for the handling of chemicals.

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Safety precautions

ISMATEC® tubing pumps are designed for pumping and dispensing applications in laboratories and industry.

As such it is assumed that Good Laboratory Practice (GLP) and our following recommendations will be observed.

- The pump must not be operated outside the destined operating and environmental conditions. In particular, it is in the user's responsibility not to exceed the specified maximum differential pressure for the pump head.
- The pump must not be used:
 - for medical applications on human beings
 - in explosion proof chambers or in the presence of flammable gases or fumes.
- The pump must be switched off when pump heads, cassettes or tubing are inserted or changed.
- The permeability of tubing depends on the material used and pressure conditions. Tubing can also become electro-statically charged. Please be aware of possible hazards when laying tubing in explosion-proof chambers.
- Pump heads consist of re-volving parts. Therefore, the pump must not be operated before the cassettes are fully snapped-in.

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Safety precautions

- Do not manipulate the pump head before the pump is switched off and disconnected from the mains supply.
- Be particularly cautious that no parts of your body such as fingers, long hair, etc. or jewellery, or loose objects such as cables or tubing, etc. can be trapped by the revolving rotor.
- Tubing can tear and burst during operation. If this could cause damage, the necessary safety measures based on the specific situation must be taken. When pumping corrosive media we recommend to place the pump in a collecting basin.

⚠ The pump should not come in contact with the liquid.

- Only new fuses, according to the specifications stated on page 8 in this manual, must be used.
- The fuse holder must not be short-circuited.
- Do not open or remove the housing while the pump is operating.
- Repairs may only be carried out by a skilled person who is aware of the hazard involved.
- For service and repairs carried out by the customer or by third-party companies ISMATEC SA denies any responsibility.

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Warranty

For all parts manufactured by
ISMATEC® 2 years

All other parts,
excluding consumables: 1 year

In case of any queries, please con-
tact your ISMATEC® representative.

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Warranty terms

We warrant the perfect func-
tioning of our products, provided
they have been installed and ope-
rated correctly according to our
operating instructions.

If production or material faults
can be proved, the defective parts
will be repaired or replaced free of
charge at our discretion.

A defective pump must be re-
turned in the original ISMATEC®
packing or in a packet of equal
quality.

The duration of the warranty is
not affected by making a claim
for warranty service.
Further claims are excluded.
Shipping costs are charged to the
customer.

Our warranty becomes invalid in the case of:

- improper operation by the
user, or if the pump is diverted
from its proper use
- unauthorized modification or
misuse by the user or by a third
party
- improper site preparation and
maintenance
- operation outside of the envi-
ronmental and electrical speci-
fications for the product
- use of third-party software,
hardware, accessories or con-
sumables purchased by the
user and which do not comply
with our specifications.

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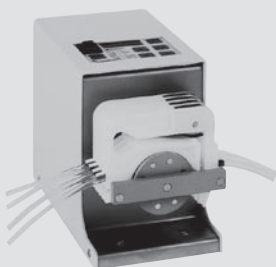
Contents of the package

- Pump REGLO *Digital* (type as ordered)
- 2 or 4 cassettes
- 1 power cord with country specific mains plug
- Operating manual

Please check the package

and its contents for transport damage. If you find any signs of damage, please contact your local ISMATEC® representative immediately.

Complaints can only be accepted within 8 days from receipt of the goods.

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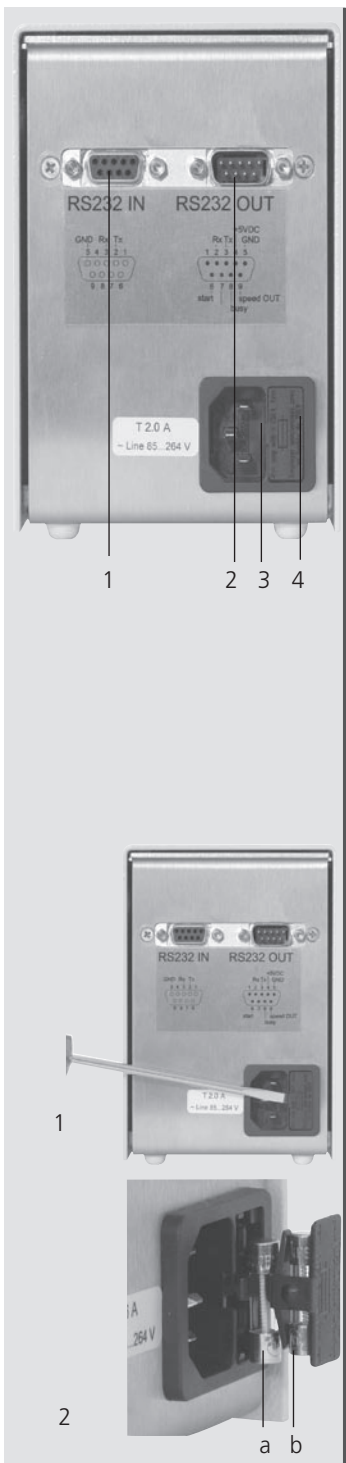
Rear panel

- 1 RS232 IN (female) for operation controlling via PC
- 2 RS232 OUT (male) for cascading up to 8 pumps
Digital input (TTL level) Run/Stop
Digital output 0 – 10 kHz (speed output), Busy (TTL Level)
- 3 Mains socket
- 4 Fuse holder

⚡ Mains voltage

Mains voltage	Fuse rating
100 – 230 V _{AC}	1 x 2.0 A 250 V slow-blow

⚠ Socket/Power cord

Use exclusively the originally supplied power cord.
The socket must be earthed (protective conductor contact).

Changing the fuses

- ➔ Switch the pump off, pull out the mains plug.
- 1 Pull out the fuse holder by opening it with a small screwdriver (size 0).
 - 2 Remove fuse (a) and replace it with the spare fuse (b) 1 x 2.0 A slow-blow.

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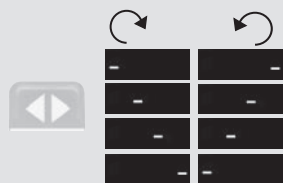
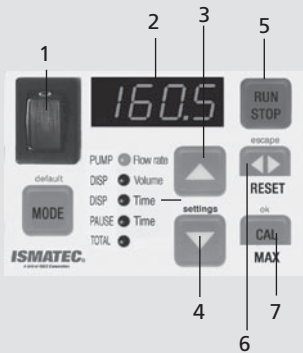
* ok

Operating panel (* = Settings page 11)

- 1 On/Off switch
- 2 Digital LED display
- 3 Increase value
- 4 Reduce value
* settings
- 5 **RUN/STOP** (double function)
 - a Starts and stops the pump
 - b Interrupts and continues a dispensing cycle (page 20)
→ The function Start/Stop can also be triggered by an external device such as a foot switch. (accessories on page 41).
- 6 **RESET** (double function)
 - a **Rotation direction**
Changes the rotation direction. By changing the rotation direction or starting the pump, the rotation direction is indicated in the display with a running hyphen.
 - b **RESET**
 - Resets the accumulated volume in mode TOTAL to zero (page 18)
 - Interrupts dispensing cycle for good (page 20)
- 7 **CAL/MAX** (double function)
 - a **CAL** = Calibrating key for the flow rate or dispensing volume
 - b **MAX** = Max. speed in the modes PUMP•Flow rate and TOTAL (keep the MAX key pressed down whilst the pump is running)

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Operating panel (* = Settings page 11)

- 8 MODE key
Selecting operating mode
* default

9 LEDs for active operating mode

- **PUMP • Flow rate**
Pumping by flow rate (ml/min)
- **DISP • Volume**
Dispensing by volume (ml)
- **DISP • Time**
Dispensing by time
- **PAUSE • Time**
Pause time for dispensing at intervals
- **TOTAL**
Pumping with Read-out of total-ly delivered volume

Start-up information

After switching on the power supply switch, the following values are displayed:

- 1 LED test »8.8.8.8.«
- 2 Firmware version
- 3 Displays the operating mode last used

⚠ **Before starting the pump for the first time or after connecting a new pump head,** the ID-code of the pump head must be entered (see basic settings, page 15).

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⚠ Please note

The settings menu requires a break of **5 seconds** between power off and on.

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Selecting the basic settings

- 1 Switch the pump off and wait for 5 seconds
- 2 Keep the settings key pressed and switch the pump on
- 3 Select the basic settings by using the ▲▼ keys
 - Cycles (dispensing cycles)
 - Back-steps (number of back-steps for volume dispensing)
 - Address (RS232, pump address)
 - Foot switch (Toggle/Direct)
 - Autostart (Standby/RUN)
 - Digits (Number of digits displayed for flow rate and dispensing volume)
 - Display Speed (Type of speed display)
 - Tube (Inner diameter)
 - Head = ID-Code (pump head)
 - Init (resets basic settings to the default values except pump head and tube diameter)
- 4 ok = confirm your selection
- 5 Select the value by using the ▲▼ keys, or
 - press escape for cancelling wrong input
 - the values can be reset to the default value by using the default key

6 ok = confirm entered value

7 escape
Quit menu basic settings

➔ **Basic settings**
refer to the following pages

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Number of dispensing steps
0..9999 (0 = Default value,
infinity)

Backsteps

Number of roller back-steps (for
drip-free dispensing in the mode
DISP Volume)
from 0 .. 100 (0 = Default value)

Address (RS232)

Pump address 1 – 8 (page 32)
(1=default value)

Foot switch

- Toggle (default value)
 - Step on foot switch
= pump starts running
 - Press foot switch again
= pump stops
- Direct
the pump runs as long as the
foot switch is pressed.

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Basic settings

Configuring Autostart

Autostart determines the pump action when the mains supply is on again, i.e. after power black-out or when the pump is controlled via mains power.

Standby

Pump is stand-by when switched on (default value)

RUN

Pumpe runs when switched on.

Digits (display)

Number of significant digits, 3 or 4 (3 = default value) when displaying flow rate and dispensing volume.

Display Speed

Display of the pump speed in mode Pump•Flow rate as

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Tubing inner diameter

Configuring the tubing inner
diameter
(0.13 mm = default)

Head Identification code (ID)
of the pump head,
see also table on page 39
(2-06 = default value)

➔ The pump head ID is com-
posed of:

- number of channels (2 or 4)
- number of rollers (6, 8 or 12)

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Basic settings

Pressing the ok key resets the following basic settings to the default values (except pump head and tube diameter):

Mode	PUMP•Flow rate
Speed	40 rpm
Dispensing volume	30 roller steps
Dispensing time	4.50 sec
Pause time	2.00 sec
Rotation direction	clockwise
Back-steps	0
Number of dispensing cycles	0
Pump head digits display	3
Foot switch	Toggle
Autostart	Standby
Display	Flow rate

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Inserting the tubing

⚠ First enter the tube i.d. as described on page 15.

- 1 ➔ Switch the pump off
- 2 Remove the cassette by slightly pressing the fixing-tongue and lifting it simultaneously
- 3 Insert the 3-stop color-coded tubing with one stopper into the cassette
- 4 Let the tubing hang down (prevent it from being twisted)
- 5 Insert the tubing with the second stopper at the other end of the cassette
- 6 Reinsert the cassette into the roller-head
- 7 ➔ Switch the pump on

➔ When the pump is idle, release all cassettes at the fixing-tongue (2).

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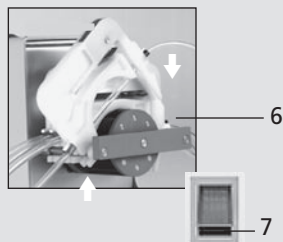
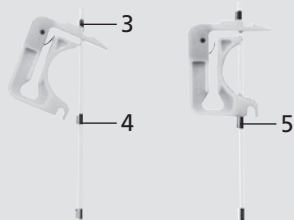
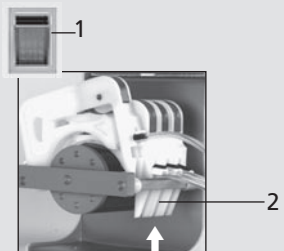
3-stop color-coded tubing

Please refer to our detailed tubing documentation.

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Pumping by flow rate

- 1 Change mode to PUMP • Flow rate
- 2 Enter the required flow rate with the ▲▼ keys (displayed in ml/min) Start pump with RUN/STOP
➔ The flow rate can also be adjusted while the pump is running.
➔ For an accurate flow rate we recommend to calibrate the pump. (see page 19).
- 3 For fast filling or emptying the system, keep pressing the CAL/MAX key while the pump is running.

➔ During the pumping process the following values can be retrieved with the MODE key:

- Flow rate (PUMP rpm)
- Totally pumped volume (TOTAL)

TOTAL volume

Change with the MODE key to TOTAL. The totally delivered volume is displayed.

Four display readings are available:

Microliter:	e.g.	17.0 µl = 17.0µ
Milliliter:	e.g.	17.0 ml = 17.0
Litre:	e.g.	17.0 l = 17.0L
Litre:	e.g.	1700 l = 1700.

For zero-setting, either

- press the RESET key once, or
- switch the pump off.

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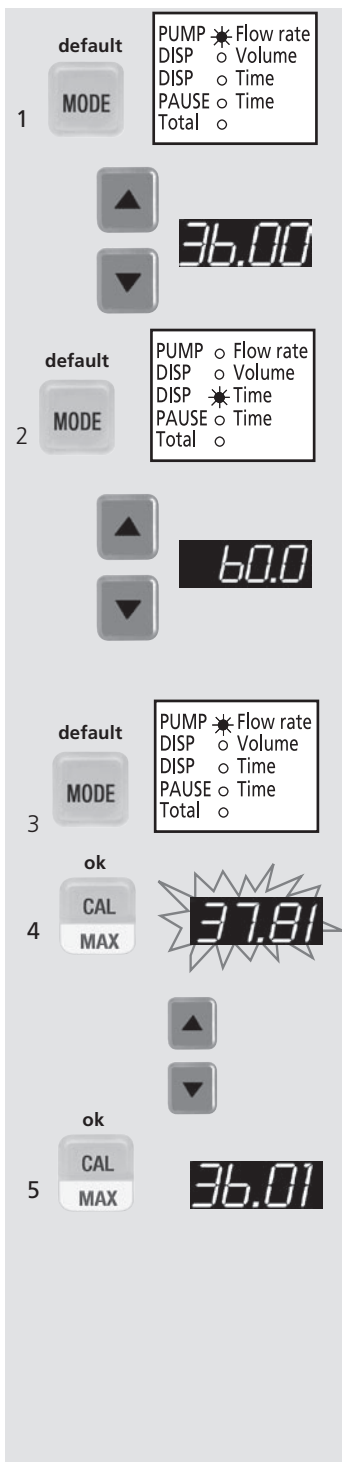
Calibrating the flow rate

- 1 Change mode to PUMP • Flow rate
Enter the required flow rate with the ▲▼ keys
➔ Check whether the correct tubing i.d. is selected (see page 15).
- 2 Change mode to DISP • Time
– Enter 60 seconds by using the ▲▼ keys
– Start pump with RUN/STOP
➔ Pump stops automatically after 60 seconds
➔ Measure the dispensed liquid by volume or weight.

Calibrate the ascertained value as follows:

- 3 Return to mode PUMP • Flow rate.
- 4 Press the CAL/MAX-button (displayed value blinks)
Enter the measured value (balance) with the ▲▼ keys
- 5 Confirm with the CAL/MAX-key (the flow rate setting is re-adjusted to the initially entered set point as close as possible). The rotation speed has a resolution of $\frac{1}{100} \text{ rpm}$. After calibration, the initially set flow rate can only be approached.

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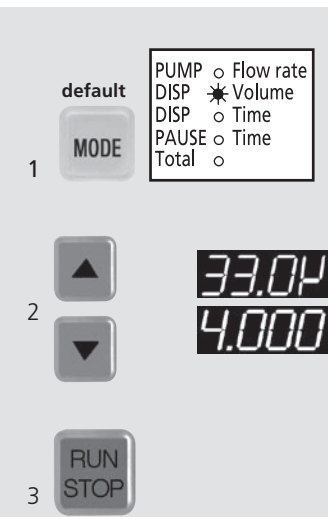
Please note

- A dispensing cycle can be interrupted by pressing the RUN/STOP key
- the volume dispensed up to now is blinking
 - by pressing the RUN/STOP key again, the dispensing cycle is continued
 - by pressing the RESET-key the dispensing cycle is stopped for good

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Dispensing by volume

- 1 Change the mode to DISP • Volume
- 2 Use the ▲▼ keys for entering the required dispensing volume. Two display readings are available:
Microlitre: e.g. 12.5 µl = **12.5µ**
Millilitre: e.g. 230.5 ml = **230.5**
➔ For dispensing accurately we recommend you to calibrate the pump (see page 21).

- 3 Start pump with RUN/STOP. The pump reduces the rotation speed shortly before the end of the dispensing cycle providing controllable and drop-precise dispensing cycles.

➔ The dispensing speed can be adjusted in the mode PUMP • Flow rate.

➔ With the ▲▼ keys the dispensing volume can also be changed **even during the dispensing process**. The newly entered set point appears briefly in the display.

➔ During the dispensing process the flow rate (PUMP • Flow rate) and the total volume can be retrieved by pressing the MODE key.

➔ The flow rate can be changed even during the dispensing process by using the ▲▼ keys.

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Calibrating the volume

- 1 Change the MODE to DISP • Volume
➔ Use the ▲▼ keys for entering the required dispensing volume
 - 2 Start with RUN/STOP
➔ The pump reduces the rotation speed shortly before the end of the dispensing cycle providing controllable and drop-precise dispensing volumes
➔ Measure the dispensed liquid by volume or weight
 - 3 Stay in the same mode (DISP • Volume) and press the CAL/MAX-key (display blinks)
 - 4 Enter the effectively dispensed volume with the ▲▼ keys
 - 5 Save with the CAL/MAX key.
➔ The calibrated volume is adjusted to the initially entered set point as close as possible. (The volume per roller step determines the smallest calibration unit. A dispensing volume is always the multiple of this unit. The initially set volume can often only be approached).
 - 6 Start pump with RUN/STOP
➔ Depending on the application conditions repeated calibration may be necessary.
- ➔ The smallest possible dispensing volume corresponds to the volume of one roller step.

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Please note

A dispensing cycle can be interrupted by pressing the RUN/STOP key

- the volume dispensed up to now is blinking
- by pressing the RUN/STOP key again, the dispensing cycle is continued
- by pressing the RESET-key the dispensing cycle is stopped for good

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Dispensing by time

The dispensing time can be entered from 0.1s – 999h.

- 1 Change the mode to DISP•Time
- 2 Enter the required dispensing time with the ▲▼ keys
 ➔ The display accelerates when pressure on the ▲ or ▼ key is maintained.

Seconds: 0.1" – 899.9" (in 0.1 s steps)
 Minutes: 15' – 899' (in 1 min steps)
 Hours: 15 h – 999 h (in 1 h steps)

- 3 Start pump with RUN/STOP

- ➔ With the ▲▼ keys the dispensing time can be changed while the dispensing cycle is in process.
- ➔ During the dispensing process the flow rate (PUMP • Flow rate) and the total volume can be retrieved by pressing the MODE-key. The flow rate can be changed during the dispensing process by using the ▲▼ keys.
- ➔ Very short dispensing times can result in non-repeatable dispensing volumes.

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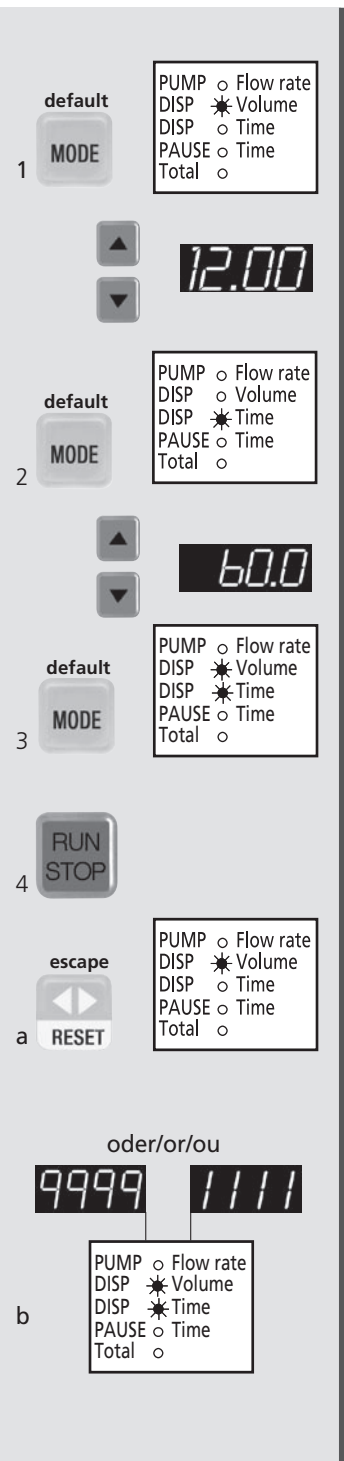
Dispensing by volume within a pre-set time

- 1 Change to mode DISP • Volume. Enter the required volume (in µl or ml) with the ▲▼ keys.
- 2 Change to the mode DISP•Time.
Enter the required time with the ▲▼ keys (0.1 s – 999 h)
- 3 Push the MODE-key once
→ Both LEDs light up simultaneously
→ 1111 and 9999 indicate an error message (see below)
- 4 Start pump with RUN/STOP
→ The pump dispenses the required volume within the pre-set time.

a Changes / Corrections

- For changing the dispensing time or dispensing volume press the RESET-key when the pump has stopped.
- The pump returns directly to MODE DISP • Volume
 - Repeat steps 1 to 3

b Error message

LEDs are blinking alternately

- 9999 = volume too large
time too short
1111 = volume too small
time too long

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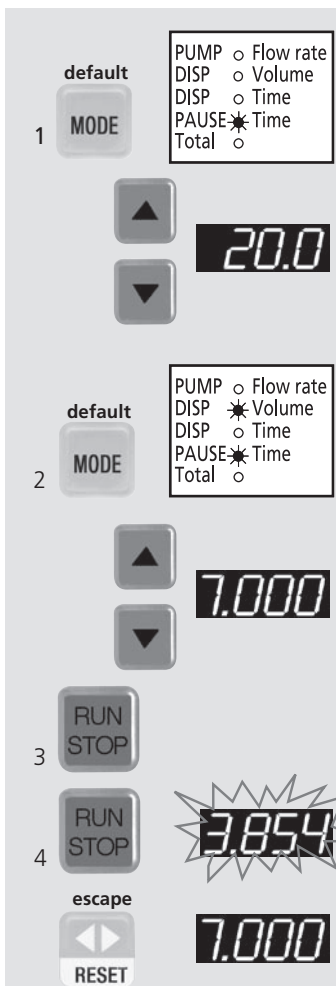
Please note

For entering the number of dispensing cycles see page 13.

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Intermittent dispensing (by volume)

Intermittent dispensing by volume with a pre-set pause time

- 1 Go to the mode PAUSE • Time.
Enter the pause time (between 0.1 s – 999 h) with ▲▼ keys.
- 2 Press MODE key until both LEDs for PAUSE and DISP • Volume light up simultaneously.
➔ Enter the dispensing volume with the ▲▼ keys.
- 3 Start with the RUN/STOP key.
➔ The pump dispenses now at intervals.
➔ Active mode blinks.
- 4 Pushing the RUN/STOP key again interrupts the dispensing process or pause period.
➔ The display blinks.
The user can now proceed as follows:
➔ By pressing the RUN/STOP key again, the dispensing or pause period is continued.
➔ By pressing the RESET-key the pump is stopped for good.

➔ With the ▲▼ keys the dispensing volume or pause time can be changed even during the dispensing process. The newly entered set point appears shortly on the display.
➔ During the dispensing process the flow rate can be retrieved with the MODE key and adjusted with the ▲▼ keys. Furthermore, the total volume can be displayed this way.

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Intermittent dispensing (by time)

Intermittent dispensing by time with a pre-set pause time

- 1 Go to mode PAUSE • Time
Enter the pause time (between 0.1 s – 999 h) with ▲▼ keys.
- 2 Press MODE key until both LEDs for PAUSE and DISP•Time light up simultaneously.
➔ Enter the dispensing time (0.1 s – 999 h) with ▲▼ keys.
- 3 Start with the RUN/STOP key.
➔ The pump dispenses now at intervals.
➔ Active mode blinks.
- 4 Pushing the RUN/STOP key again interrupts the dispensing process or pause period.
➔ The display blinks.

The user can now proceed as follows:

- ➔ By pressing the RUN/STOP key again, the dispensing or pause period is continued.
- ➔ By pressing the RESET-key the pump is stopped for good.

- ➔ With the ▲▼ keys the dispensing or pause time can be changed even during the dispensing process. The newly entered set point appears shortly on the display.
- ➔ During the dispensing process the flow rate can be retrieved with the MODE key and adjusted with the ▲▼ keys.

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Please note

With RUN/STOP a dispensing cycle can be interrupted and resumed.

With RESET an interrupted dispensing cycle can be definitely stopped.

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Number of dispensing cycles

The number of dispensing cycles can be entered when dispensing at intervals (by time and volume).

See basic settings on page 13.

➔ During the pauses the display shows the remaining number of dispensing cycles.

➔ During long pauses, the remaining interval time can be checked with the CAL/MAX key and adjusted by using the ▲▼ keys.

➔ At any time, the dispensing process can be interrupted by pushing the RUN/STOP-key. When resuming the process, the pump continues with the remaining number of dispensing cycles.

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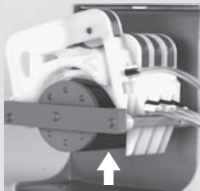
Pumping against pressure

The REGLO *Digital* can be used for continuous duty at a differential pressure of max. 1.0 bar (smaller tube sizes and/or lever type cassettes allow higher pressures).

When the pump is idle

When the pump is idle, we recommend you to release the tubing from pressure. Releasing the cassette on the right side is sufficient. This helps to protect the tubing from unnecessary strain and prolongs its service-life.

⚠ Syphoning effect

When the tubing is released from squeezing, the fluid can flow back to the reservoir.

Overcurrent protector

The drive REGLO *Digital* features an overload protector. It disables the pump drive in case of a motor blockage. When an overload condition occurs, it is indicated in the display by the letters 'OL' and the pump is stopped.

Whenever this situation occurs, **the pump must be switched off immediately**. Wait for 5 sec. before switching on the pump again. Before the pump is re-started, it is most important to check the reason for the overload (eg. too high differential pressure etc.).

➔ Only when the cause of the overload has been detected and the failure corrected accordingly may the pump be started again.

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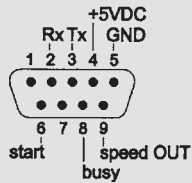
Please note

The analog interface is integrated in the port 'RS232 OUT' and can be used for remote control (e.g. by foot-switch).

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Analog interface

Pin 4: +5VDC

For connecting an external control unit (current <100mA).

Pin 5: GND

Point of reference for all signals described.

Pin 6: start

Connection for an external trigger device (foot-switch, etc.) For operating modes, refer to Basic settings "Foot switch" page 13

Pin 8: busy

Informs about pump status:

Pumpe in operation
= output +5V
Pumpe in stand-by
= output on ground
(The output resistance is 4.7 kΩ)

Pin 9: speed OUT

A square-wave signal with a frequency proportional to the speed (rpm) is at your disposal:

100 rpm	5000 Hz
160 rpm	8000 Hz

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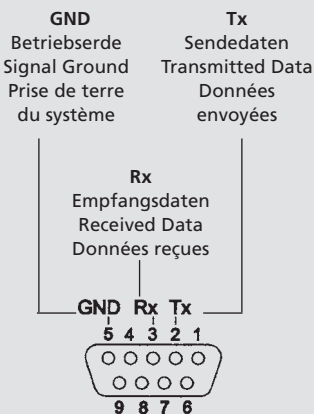
Please note

When switching on the REGLO *Digital* pump, the serial interface will only respond after a delay of 5 seconds.

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Serial interface

RS232 IN (female)

A 9-pin D-socket is available on the rear panel of the pump.

Setting the pump address

Each command string begins with the pump address (1 – 8).

It is factory-set at 1 and can be changed with @.

The address will be stored permanently (even after the pump has been switched off).

The assignment of the address enables the user to control up to 8 ISMATEC® drives via one interface.

Structure of the commands

The address is followed by a character. Some commands have an additional parameter which always consists of 4 or 5 figures.

The command string **is completed by** a carriage return (ASCII 13). The pump confirms most of the commands with an asterisk *. Yes/No inquiries are answered by + (yes) or – (no).

Multi-digit replies are concluded by (cr) ASCII 13 and (lf) ASCII 10.

Incorrect command strings are answered by #.

Numerical values are confirmed as 3 to 5-digit figures. Four of the five digits are numerals, one digit is either a decimal point or a preceding blank space.

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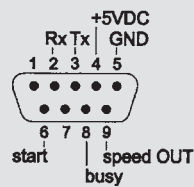
RS232 OUT (male)

Pin 2: RS232 Rx

Is used for connecting additional pumps together with pin 3/5.

Pin 3: RS232 Tx

Is used for connecting additional pumps together with pin 2/5.

Pin 5: GND

Point of reference for all signals described.

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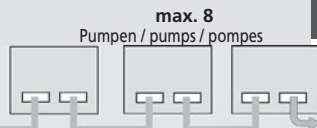
**Each pump requires
an individual address.
(see page 31)**

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Cascading several pumps

Providing that an appropriate software is available, up to 8 ISMATEC® pumps can be controlled via one PC.

In order to establish a cascade you need to connect the pumps with the connecting cable (AG0013) as shown on the left.

The same cable AG0013 serves to connect the first pump to the PC, provided that it has a 9-pin RS232 connector.

If the PC is equipped with a 25-pin connector, you need an additional adapter XC0009.

- ➔ Each pump must be allocated its own individual address.

For this procedure, you first have to connect each pump in turn to the RS232 interface of the PC. Two or more connected pumps must not have the same address. Only when all the pumps have been assigned their addresses can you cascade them as shown on the left.

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Please note

The values indicated in this table are only approximate and determined as follows:
ml/min, per channel, with water and Tygon® tubing, without differential pressure

Order No. = Tygon® ST (R3603)

We will be pleased to help the user to select the tubing material. However, the user himself has the final responsibility for the selection of the correct tubing material. For reproducible results we recommend you to use tubing from ISMATEC®.

Please ask for the detailed tubing documentation.

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
REGLO Digital/ISMATEC SA / 13.09.07/CMC

Artisan Technology Group - Quality Instrumentation ... Guaranteed | (888) 88-SOURCE | www.artisantg.com

39

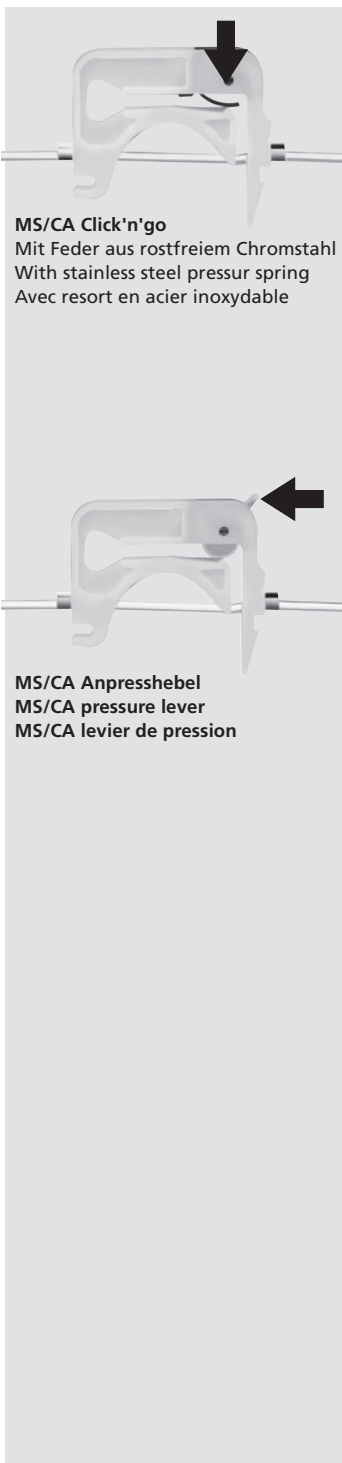MS/CA Click'n'go
Mit Feder aus rostfreiem Chromstahl
With stainless steel pressur spring
Avec ressort en acier inoxydable

MS/CA Anpresshebel
MS/CA pressure lever
MS/CA levier de pression

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Spare cassettes

Cassette	MS/CA Click'n'go
Material	POM-C
Order No.	IS 3510
Cassette	MS/CA pressure lever
Material	POM-C
Order No.	IS 0649
Cassette	MS/CA pressure lever
Material	PVDF
Order No.	IS 3629

➔ The cassettes with pressure lever are available on request. This type of cassette may provides better results both for the Tygon MH tubing and at elevated differential pressure conditions.

➔ The cassettes allow the user to insert tubes with different diameters and materials on the same roller-head.

Cassette material

Materials POM-C

Polyoxymethylene-Copolymer
– good chemical resistance to many organic solvents and strong alkaline chemicals
– is affected by strong acids and oxidizing substances
– UV-stabilized and stable up to temperatures of 80°C/176°F (dry, continuous use) or 136°C/277°F (dry, for a short time)

Materials PVDF

Polyvinylidene fluoride
– very good chemical resistance to acids and most aliphatic, aromatic and chlorinated solvents
– not suitable for long contact with esters, ketones, amines and strong alkaline chemicals
– stable both to UV radiation and temperatures up to 110°C/230°F (continuous use) or 142°C/288°F (for a short time)

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Preparing the cassettes with pressure lever

- 1 ➔ Switch the pump off
- 2 Push the pressure lever down (it must be relaxed)
Insert the tubing with 2 colour-coded collars
- 3 Reinsert the cassette(s) onto the rollers. Ensure that it snaps in on both sides.
- 4 ➔ Switch the pump on and start it
- 5 Lift the pressure lever step by step until the medium is drawn in and a continuous flow is achieved (no further step than necessary!)

Note

- Under certain circumstances (e.g. new tubing, elevated differential pressure, etc.) we recommend the user to lift the lever by just one more step
- Never over-tighten the pressure lever
- The setting of the pressure lever determines the occlusion
- The optimum occlusion depends on tubing material, inner diameter, wall thickness, age of tubing, and application conditions
- As a side-effect the occlusion setting also affects the flow rate. This, however, should not encourage the user to use the pressure lever for flow rate adjustments because long-term impact is not predictable
- If, however, the occlusion setting is used to adjust the flow rate, the channel-to-channel flow rate conformity may lead to good results in the beginning of the operation, yet the individual flow rates may increasingly diverge at a later stage
- For an optimum long-term flow rate stability and channel-to-channel conformity we suggest using the Click'n'go cassettes the flow rate

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Service and repairs

Provided the REGLO *Digital* is operated properly and in compliance with this manual, the tubing is the only part that is subject to wear and tear.

For repairs send the complete pump with detailed description of the failure to your ISMATEC® distributor. Please use the original ISMATEC® packing or a packet of equal quality.

Replacement parts

For repairs you intend to carry out by yourself (after the warranty period), you can ask from your ISMATEC® distributor:

- Parts lists
- Replacement parts for pump drive
- Wiring diagrams

Please give information on defect, date of purchase, serial No., and model.

Disposal

Please retain packing materials until the product warranty ends. Afterwards please discard packing materials in an environment friendly manner according to local regulations.

Once the useful life of the product has ended, please ensure proper disposal according to local laws. Plastic and electronic components should be disposed of at a recycling facility. Please refer to local regulations regarding proper disposal.

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⚠

Please note

We also recommend you to observe our Warranty Terms as well as our Terms and Conditions of Sale.

In case of any queries, please contact your local ISMATEC® representative.

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2 channels1.6–160 rpm

4 channels1.0–100 rpm

digitally adjustable in steps of 0.01 rpm

Differential pressure

max. 1.0 bar (14.5 psi)

Remote control

via RS232 and analog interface (see pages 30–38)

Mains connection

100 – 230 V_{AC} / 50 – 60 Hz

Fuse rating power supply

2.0 A / 250 V, slow-blow

Power consumption

max. 20 W

Protection rating

IP 30

Operating conditions

Temperature+41 to +104°F

Rel. humiditymax. 80%

– not condensing, at normal laboratory conditions

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
Ask for the latest
ISMATEC catalogue.

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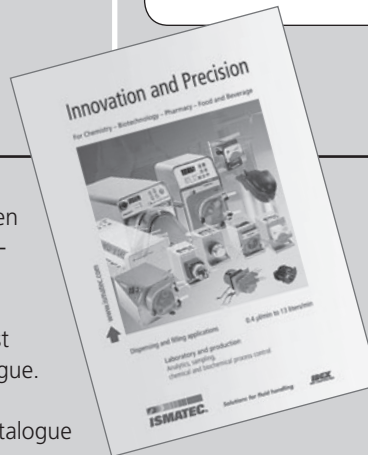ISMATEC® - Your competent partner for demanding metering and dispensing applications.

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Laboratoriumstechnik GmbH