

TERUMOBCT

Unlocking the Potential of Blood

TSCD® II STERILE TUBING WELDER

THE ABSOLUTE REFERENCE IN STERILE CONNECTION



SMOOTH, FAST AND EFFICIENT

TSCD II STERILE TUBING WELDER

Information Panel

Backlight LCD provides status and operational prompts.



With two decades of experience and extensive use, Terumo BCT's sterile connecting devices have become indispensable for accomplishing quick and safe tubing connections that maintain a functionally closed system. This level of quality and reliability is assured, regardless of wet or dry tubing connections. The new TSCD provides a new level of speed and assurance.

TSCD II: Strong, Consistent Welds

All Terumo BCT's systems are subject to a minimum tensile strength requirement.

TSCD Wafers

The TSCD II uses the same wafers as the TSCD. One box includes two cassettes; each cassette contains 70 single-use wafers. Cross-contact or contamination is avoided because each wafer is discarded after a single use. (Product code: SC*W017)



Automatic Wafer Feed

Wafers are automatically replaced after each use. When the wafer cassette empties or the wafer disposal box is full, the TSCD II displays a warning message.

Clamps

The clamping system is designed to include wet to wet applications. A robust clamp mechanism makes TSCD II the first choice for intensive use. The interconnected left and right clamping mechanism facilitates connections and helps to ensure that the TSCD II is the primary choice for efficiency and ease of use.

Wafer Disposal Box

For safety and convenience, the used wafers eject into a built-in disposal box. The disposal box holds as many as 300 used wafers, and an integrated sensor monitors its volume at all times. These features help prevent wafer disposal box overflow.

Portable

TSCD II is designed for mobility. Its lightweight, compact size and pop-up handle make it easy to take with you.



Data processing system available.

APPLICATIONS INCLUDE:


- Component Pooling
- Leukoreduction
- Collection Set Modification
- Component Aliquoting
- Apheresis Set Modification
- Quality Control Sampling
- Stem Cell Processing
- Cell Washing and Freezing
- Other Applications Involving Closed System PVC Tubing Connections




TSCD II Terumo BCT Sterile Tubing Welder

Tube Cutting/Welding Process 300°C (572°F) Wafer Temperature


Terumo BCT's patented wafers heat up to 300°C (572°F) – a temperature high enough to maintain sterility during cutting and welding without allowing particulates or chemical residues to form in the tubes.

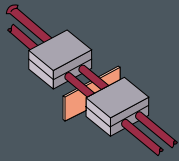


The TSCD wafer is heated to approximately 300°C (572°F)

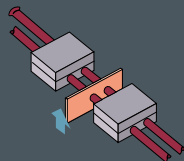

The wafer moves up and melts through both pieces of tubing

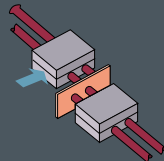

The left clamp moves horizontally until the tubings are aligned

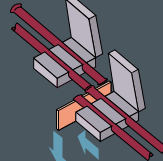

When the wafer is retracted, the right clamp pushes the tubing together

Terumo BCT Service Center

To ensure system reliability, Terumo BCT's Service Center conducts a series of strict, in-house operational tests before shipment.

Specifications

Name of Product	Terumo BCT Sterile Tubing Welder TSCD®-II
Product Code	ME-SC 203A
Tubing Material	Polyvinyl chloride (PVC)
Tubing Sizes	O.D.: 3.9 – 4.5 mm / I.D.: 2.9 – 3.1 mm
Wafer	TSCD® wafer / Product code: SC*W017
Operating Conditions	Ambient temperature: 10°C to 40°C (50°F to 104°F) Relative humidity: 10% to 80%, non-condensing Maximum altitude: 2,000 m (6,562 ft)
Storage Conditions	Ambient temperature: -10°C to 50°C (14°F to 122°F) Relative humidity: 10% to 95%, non-condensing
Power Supply	AC 100V to 240V±10%, 50/60 Hz
External Dimensions	224 (W) x 177 (H) x 342 (D) mm
Weight	Approx. 6.5 kg (14.33 lbs)
Accessories with TSCD II	AC power cable, bag supports, fuses (250V/4.0A), cassette of TSCD® wafers, replacement air filter and operating instructions
Communication Ports	RS-232 for TSCD® IS (TSCD® Information System) and barcode reader, 2 RS-485 for TSCD® IS (for E.U. only)



Tensile Strength Test



Pressure Test

* Tubing shall meet the specifications for the welder. Tubing outside of the specifications and combinations of different tubing shall be subject to validations. Please contact your Terumo BCT representative for details.

Manufactured by: Terumo Europe N.V., 3001 Leuven, Belgium
Distributed by: Terumo BCT

TERUMOBCT

Terumo BCT, Inc.

10811 West Collins Ave.
Lakewood, Colorado 80215-4440
USA
USA Phone: 1.877.339.4228
Phone: +1.303.231.4357
Fax: +1.303.542.5215

Terumo BCT Europe N.V.

Europe, Middle East and Africa
Ikaroslaan 41
1930 Zaventem
Belgium
Phone: +32.2.715.05.90
Fax: +32.2.721.07.70

Terumo BCT (Asia Pacific) Ltd.

Room 3903-3903A, 39/F
ACE Tower, Windsor House
311 Gloucester Road
Causeway Bay, Hong Kong
Phone: +852.2283.0700
Fax: +852.2576.1311

Terumo BCT Latin America S.A.

La Pampa 1517 - 12th Floor
C1428D2E
Buenos Aires
Argentina
Phone: +54.11.5530.5200
Fax: +54.11.5530.5201

Terumo BCT Japan, Inc.

Ebisu South One Bldg, 9F
1-7-8, Ebisu-minami,
Shibuya-ku, Tokyo 150-0022
Japan
Phone: +81.3.6743.7890
Fax: +81.3.6743.9800