

Specification Sheet

Liquid Chromatograph Mass Spectrometer

LCMS-8040

Triple Quadrupole Mass Spectrometry is the method of choice for definitive identification and reproducible quantification of trace level analytes in complex samples for a variety of applications such as pharmacokinetics in pharma/CRO, food, environmental, forensic/toxicology and newborn screening. Combined with chromatographic resolving power of our world leading UHPLC systems and maintaining Shimadzu's proprietary ultrafast technologies (UF Technologies) which includes high-speed MRM transitions, MS/MS acquisitions, and ultra-high speed polarity switching, the LCMS-8040 can dramatically improve analytical throughput with ultra-high speed performance. In addition, the newly improved ion optical system, UF-Lens and unique UFsweeper II collision cell technology yield higher sensitivity and expand the potential range of LC/MS/MS applications.



Instrument

Model	LCMS-8040
Mass range	m/z 2 to 2,000
Sensitivity	ESI positive: 1 pg reserpine, S/N>30,000:1 (RMS) ESI negative: 1 pg chloramphenicol, S/N>30,000:1 (RMS)
Resolution	$R < 0.7 \mu$ FWHM
Mass stability	0.05 μ / 12 hr
Cross talk	< 0.003%
Minimum pause time	1 msec
Minimum dwell time	0.8 msec
Scan speed	Max 15,000 μ /sec (0.1 μ step: 150,000 data points/sec)
Polarity switching time	15 msec
Interface	ESI (Standard), APCI (Optional), DUIS (Optional)
Applicable LC flow rate	ESI 1 μ L/min to 2 mL/min
MRM transition speed	Max 555 channels/sec
DL maximum temperature	300 °C
Block heater maximum temperature	ESI/DUIS: 500 °C APCI: 300 °C
Interface maximum temperature	APCI: 500°C

Analysis mode	Q1 Scan/SIM Q3 Scan/SIM MRM Precursor ion scan Product ion scan Neutral loss scan
---------------	--


Mass Analyzers and Detector

Mass analyzers	Q1 & Q3 are molybdenum hyperbolic mass filters with pre-rods; Q1 includes post-rods
Collision cell	Tapered multipole type ultra-high-speed collision cell (UFSweeper II collision cell)
Detector	Secondary electron multiplier with off-axis conversion dynode
Ion optics	Q-array focus optics operating in Field-Flow mode, multipole transfer optics
Digital detection system	Operates in pulse counting mode for fastest operation
Detection mode	Ultra-fast positive/negative ion switching
Dynamic range (pulse counting)	2×10^7 cps
Vacuum system	Rotary pump: 1 unit Vacuum pumping speed: 28 m ³ /hr Triple-inlet turbo molecular pump: 1 unit 40 L/sec, 260 L/sec, 210 L/sec

Software

Workstation	LabSolutions LCMS
Operation system	Windows 7 32 Bit, 64 Bit
Instrument control	Prominence and Nexera series
MS acquisition mode	Scan (Max. 1,000 events), SIM (Max. 1,000 events × 32 channels)
MS/MS acquisition mode	MRM (Max. 1,000 events × 32 channels) Product ion scan Precursor ion scan Neutral loss scan
Auto-tuning	Possible to optimize sensitivity and resolution in both positive and negative ionization mode as well as mass calibration

Installation Example

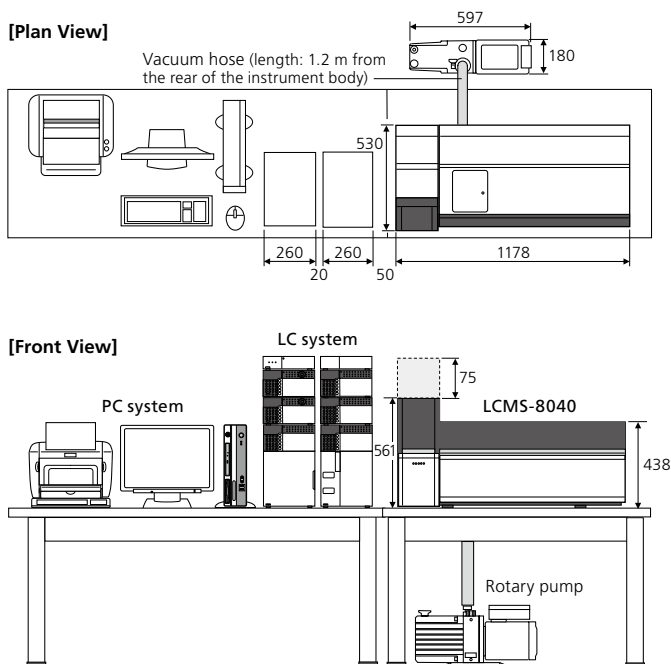

Units: mm

Installation Conditions

Temperature	18 to 28 °C
Humidity	20 to 70 % (Non-condensing)
Size	1,178 mm (W) × 530 mm (D) × 636 mm (H)
Weight	130 kg
Power supply	MS unit: AC 230 V 15 A (50/60 Hz) Single-phase
Gas requirements	Nitrogen gas: Maximum consumption 25 L/min, Purity greater than 97 % Argon: Purity greater than 99.99 % as CID gas

The above are not standard installation specification. All LCMS-8040 instruments will be installed and tested in accordance with standard performance tests as detailed in the Shimadzu document ZEAH-0450, Shimadzu High-Performance Liquid Chromatograph Mass Spectrometer LCMS-8040 Installation Standard.

This product is for research use only, not for use in diagnostic procedures.



Shimadzu Corporation

www.shimadzu.com/an/

For Research Use Only. Not for use in diagnostic procedures.

This publication may contain references to products that are not available in your country. Please contact us to check the availability of these products in your country.

Company names, products/service names and logos used in this publication are trademarks and trade names of Shimadzu Corporation, its subsidiaries or its affiliates, whether or not they are used with trademark symbol "TM" or "®". Third-party trademarks and trade names may be used in this publication to refer to either the entities or their products/services, whether or not they are used with trademark symbol "TM" or "®". Shimadzu disclaims any proprietary interest in trademarks and trade names other than its own.

Shimadzu does not assume any responsibility or liability for any damage, whether direct or indirect, relating to the use of this publication.

© Shimadzu Corporation, 2017